


NEONATOLOGY

ENTERAL

Syringe pump for enteral nutrition

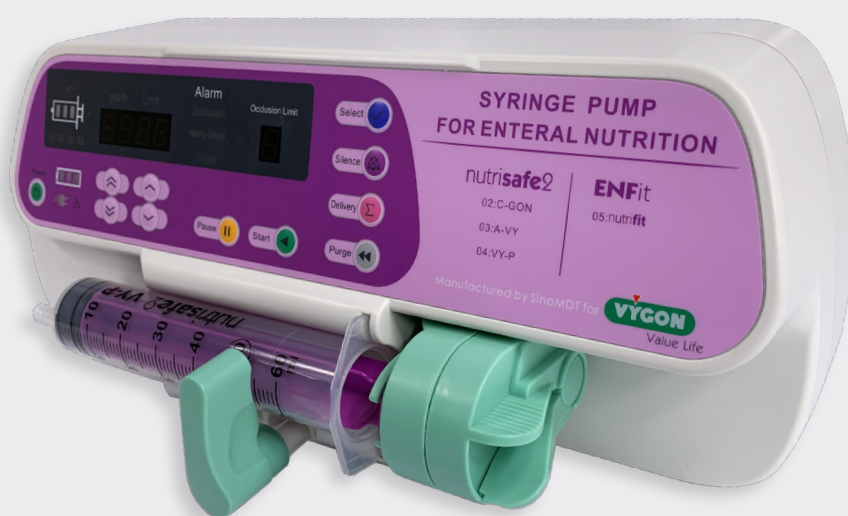

ACCURACY IN ENTERAL FEEDING

Designed to meet all patients needs in stationary medical care, this easy-to-use high-quality syringe pump is dedicated to accurate enteral feeding.



The following syringes are already calibrated on the device. Others channels are available to calibrate new syringe brands (contact your biomedical department).

| Channel Syringe vol. | nutrisafe2 | | | nutrifit |
|----------------------------|-------------------|------------------|------------------|----------------------|
| | Channel 2 / C-GON | Channel 3 / A-VY | Channel 4 / VY-P | Channel 5 / nutrifit |
| 10ml | 1015.102 | 1015.103 | 1015.104 | 1015.107 |
| 20ml | 1015.202 | 1015.213 | 1015.204 | 1015.207 |
| 30 / 35ml | 1015.352 | - | - | - |
| 50ml | 1015.602 | 1015.603 | 1015.604 | 1015.607 |

PUMP FUNCTIONS

- **Automatic syringe size identification**
- **Optic and sound alarms**
- **Possibility to adjust sound level**
- **Up to 2000 history records**
- **Delivered volume inquiry**
- **Occlusion detection pressure setting**
- **Alarms:**
 - Occlusion
 - Prior administration end
 - Administration end
 - Battery low
 - Out of battery
 - Mains failure
 - Equipment error
 - Flowrate overlimit
 - Administration completion
 - No operation within 2 minutes

ADVANTAGES



Compatible

For all brands of nutrisafe2 and nutrifit syringes



Powerful

7 hour battery autonomy



Accurate

Accuracy of flowrate including syringe $\pm 2\%$

TECHNICAL FEATURES

| | | |
|---------------------------------------|---|---|
| Dimension (LxWxH mm) | 306mm x 140mm x 146mm | |
| Weight | 2.25kg (with clamp base) | |
| Accuracy of flow rate (Nominal value) | Accuracy including syringe: within $\pm 1\%$ Mechanical accuracy: within $\pm 1\%$ | |
| Purge rate | 1500 mL/h 900 mL/h 600 mL/h 400 mL/h | 50 mL syringe 30/35 mL syringe 20 mL syringe 10 mL syringe |
| Flow rate | 0.1 to 1500 mL/h 0.1 to 900 mL/h (0.1 mL/h step) 0.1 to 600 mL/h (0.1 mL/h step) 0.1 to 400 mL/h (0.1 mL/h step) | 50 mL syringe 30/35 mL syringe 20 mL syringe 10 mL syringe |
| Increments | Specified value 0.1 to 999.9 mL/h, 0.1 mL/h step Specified value 1000 mL/h or above, 1mL/h step | |

PRODUCT CONFIGURATION

Product code: SN-50C6

- Syringe pump for enteral nutrition
- Power cords (Europe (plug C) / UK (plug G) / US (plug A))
- User manual and sticker
- Packing list
- Certificate of compliance
- Warranty card

11/22 - DB DIGE 12 163E - V3

FOR FURTHER INFORMATION, PLEASE CONTACT: questions@vygon.com

The specifications shown in this leaflet are for information only and are not, under any circumstances, of a contractual nature.

VYCON – 5, RUE ADELIN – 95440 ECOUEN - FRANCE

RECEPTION: +33 (0)1 39 92 63 63 – SERVICE CLIENTS FRANCE: +33 (0)1 39 92 63 81

EXPORT CUSTOMER SERVICE: +33 (0)1 39 92 64 15

www.vygon.com



Value Life